

Product Description

— A helping hand for your research

Active Arl3

Catalog Number: 26925 Gene Symbol: ARL3

Description: Anti-active Arl3 Mouse Monoclonal Antibody

Background: Arl3 (Arf-like 3) is an Arf family protein. Arl3 binds microtubules to alter specific aspects of cytokinesis. It has been proposed that retinitis pigmentosa 2 functions in concert with Arl3 to link the cell membrane and the cytoskeleton in photoreceptors as part of the cell signaling or vesicular transport machinery.

Immunogen: Recombinant Arl3 protein (active

form)

Tested applications: IP, IHC, IF

Recommended dilutions:

1 μg for 1~2 mg total cellular proteins for IP

1:50~100 dilutions for IHC and IF

Concentration: 1 mg/ml

Host: Mouse

Clonality: Monoclonal
Purity: Purified from ascites

Format: Liquid Storage buffer: Preservative: no

Constituents: PBS (without Mg²⁺ and Ca²⁺), pH

7.4, 150 mM NaCl, 50% glycerol.


Species Reactivity: Arl3 antibody recognizes

active Arl3-GTP from vertebrates.

Storage Conditions: Store at -20°C. Avoid

freeze / thaw cycles.

Immunoprecipitation and Western blot:

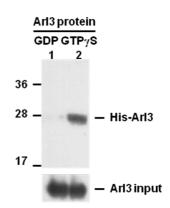

IP: anti-active Arl3 mAb
IB: anti-Arl3 mAb

Arl3 Activation Assay. Purified His-tagged Arl3 proteins (Cat. #10152) were immunoprecipitated with the anti-active Arl3 monoclonal antibody (Cat. #26925) after treated with GDP (lane 1) or GTPγS (lane 2), and was blotted with anti-Arl3 monoclonal antibody (Cat. #21070). Input controls are shown in the bottom panel.

Tel: (610) 945-2007 Fax: (610) 945-2008 Web: <u>www.neweastbio.com</u>