

Product Description

A helping hand for your research

CTNNB1 (S37Y)

Catalog Number: 26296 Gene Symbol: CTNNB1

Description: Anti-CTNNB1 (S37Y) Mouse Monoclonal Antibody

Background: CTNNB1 protein is a dual function protein. It is a subunit of a complex of proteins that from adherens junctions, which are important for the establishment and maintenance of epithelial cell layers by regulating cell growth and adhesion between adjacent cells. CTNNB1 protein also pulls double duty as an intracellular signal transducer in the Wnt signaling pathway. Mutations has been implicated in the pathogenesis of several cancers.

Immunogen: A synthetic peptide from the internal region of CTNNB1 which includes the mutation of S37Y, human origin.

Tested applications: ELISA, WB, IHC

Recommended dilutions:

ELISA: 1:1000-1:5000 WB: 1:500-1:2000 IHC: 1:50-1:100 Concentration: 1 mg/ml

Host: Mouse

Clonality: Monoclonal

Purity: Purified from ascites

Format: Liquid **Storage buffer:** Preservative: no


Constituents: PBS (without Mg²⁺ and Ca²⁺), pH

7.4, 150 mM NaCl, 50% glycerol

Species Reactivity: Recognizes S37Y mutant, but not wild-type CTNNB1 of vertebrates. **Storage Conditions:** Store at -20°C. Avoid

freeze / thaw cycles.

Western blot:

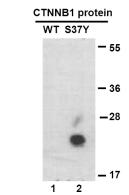

WB: anti-CTNNB1 (S37Y) mAb

Western blot analysis of recombinant CTNNB1 (S37Y) and wildtype proteins.

Purified His-tagged CTNNB1 (S37Y) protein (lane 2) and corresponding wildtype protein (lane 1) were blotted with anti-CTNNB1 (S37Y) monoclonal antibody (Cat. #26296).

Immunofluorescence:


Tel: (610) 945-2007 FAX: (610) 945-2008

Web: www.neweastbio.com



Product Description

— A helping hand for your research

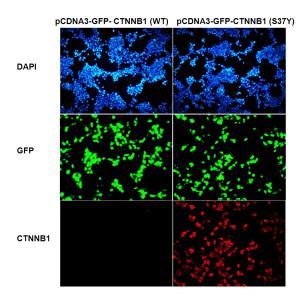

Immunofluorescence of cells expressing CTNNB1 proteins with anti-CTNNB1 (S37Y) antibody. HEK293T cells were transfected with pCDNA3-GFP-CTNNB1 (WT) plasmid (left column) or pCDNA3-GFP-CTNNB1 (S37Y) plasmid (right column), then fixed and stained with anti-CTNNB1 (S37Y) monoclonal antibody (Cat. #26296).