

## BRAF (D594G)

**Catalog Number:** 26200

**Gene Symbol:** BRAF, BRAF1, RAFB1

**Description:** Anti-BRAF (D594G) Mouse Monoclonal Antibody

**Background:** BRAF belongs to a family of serine-threonine protein kinases. As a part of a signaling pathway known as the the RAS/MAPK pathway, it plays an important part in many cellular processes. The processes include cell proliferation, differentiation and transcriptional regulation. Mutations in the BRAF gene cause diseases. Inherited mutations in BRAF cause cardiofaciocutaneous syndrome. Acquired mutations in BRAF have been found in cancers.

**Immunogen:** A synthetic peptide from the internal region of BRAF which includes the mutation of D594G, human origin.

**Tested applications:** ELISA, WB

**Recommended dilutions:**

ELISA: 1:1000-1:5000

WB: 1:500-1:1000

**Concentration:** 1 mg/ml

**Host:** Mouse

**Clonality:** Monoclonal

**Purity:** Purified from ascites

**Format:** Liquid

**Storage buffer:**


Preservative: no

Constituents: PBS (without Mg<sup>2+</sup> and Ca<sup>2+</sup>), pH 7.4, 150 mM NaCl, 50% glycerol

**Species Reactivity:** Recognizes D594G mutant, but not wild-type BRAF of vertebrates.

**Storage Conditions:** Store at -20°C. Avoid freeze / thaw cycles.

### Western blot:

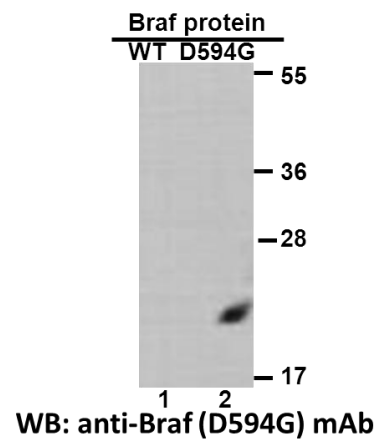

**Western blot analysis of recombinant BRAF (D594G) and wildtype proteins.** Purified His-tagged BRAF (D594G) protein (lane 2) and corresponding wildtype protein (lane 1) were blotted with anti-BRAF (D594G) monoclonal antibody (Cat. #26200).

For research use only. Not for diagnostic or therapeutic applications.



# Product Description

— *A helping hand for your research*

---

---

**For research use only. Not for diagnostic or therapeutic applications.**

**24 Whitewoods Lane, Malvern, PA 19355    Tel: (610) 945-2007    FAX: (610) 945-2008**  
**E-mail: [sale@neweastbio.com](mailto:sale@neweastbio.com)    Web: [www.neweastbio.com](http://www.neweastbio.com)**