

Product Description

A helping hand for your research

BMP3

Catalog Number: 26051 Gene Symbol: BMP3

Description: Anti-BMP3 Mouse Monoclonal Antibody

Background: BMP3 belongs to the TGF-beta family. It negatively regulates bone density and antagonizes the ability of certain osteogenic BMPs to induce osteoprogenitor differentiation and ossification. BMP3 is secreted and highly expressed in fracture tissue, particularly in osteoblasts, osteoclasts and chondroblasts.

Immunogen: A synthetic peptide near the

N-terminus of BMP3, human origin.

Tested applications: ELISA, WB, IHC, IP

Recommended dilutions:

ELISA: 1:1000-1:5000 WB: 1:1000-1:5000 IHC: 1:50-1:100 Concentration: 1 mg/ml

Host: Mouse

Clonality: Monoclonal

Purity: Purified from ascites

Format: Liquid **Storage buffer:** Preservative: no

Constituents: PBS (without Mg²⁺ and Ca²⁺), pH

7.4, 150 mM NaCl, 50% glycerol


Species Reactivity: Anti-BMP3 antibody

recognizes BMP3 of vertebrates.

Storage Conditions: Store at -20°C. Avoid

freeze / thaw cycles.

Western blot:

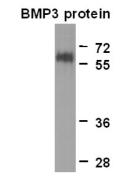

WB: Anti-BMP3 mAb

Western blot analysis of recombinant BMP3 protein. Purified His-tagged BMP3 protein (Cat. #10194) was blotted with anti-BMP3 monoclonal

Tel: (610) 945-2007 FAX: (610) 945-2008

Web: www.neweastbio.com

antibody (Cat. #26051).