

# **Product Description**

A helping hand for your research

Product Name: Rab7 Protein

Synonyms: Member RAS oncogene family, Rab7a,

MGC102153

Catalogue Number: 10167

Source: Human, recombinant full length, His6-tag

Expression Host: E. coli Molecular Weight: 23 kDa Purity: >99% by SDS-PAGE

Introduction: Members of the Rab family of Ras-related GTP-binding proteins are important regulators of vesicular transport and are located in specific intracellular compartments. Rab7 has been localized to late endosomes and shown to be important in the late endocytic pathway. In addition, it has been shown to have a fundamental role in the cellular vacuolation induced by the cytotoxin VacA of Helicobacter pylori.

# **Amino Acid Sequence** (1-207)

MTSRKKVLLKVIILGDSGVGKTSLMNQYVNKKFSNQYK
ATIGADFLTKEVMVDDRLVTMQIWDTAGQERFQSLGVA
FYRGADCCVLVFDVTAPNTFKTLDSWRDEFLIQASPRD
PENFPFVVLGNKIDLENRQVATKRAQAWCYSKNNIPYF
ETSAKEAINVEQAFQTIARNALKQETEVELYNEFPEPI
KLDKNDRAKASAESCSC

#### **Properties**

Physical Appearance (form): Dissolved in 20mM

Tris-HCl, pH8.0, 150mM NaCl.


Physical Appearance (form): White or clear

Concentration: 1 mg/mL

Storage: -80°C

## **Preparation Instructions**

Adding of 10 mM  $\beta$ -mercaptoethanol or 1 mM DTT into the solution to protect the protein is recommended and using of non-ionic detergents such as n-Dodecyl  $\beta$ -D-maltoside (DoDM) or polyethylene detergents (e.g.  $C_{12}E_{10}$ ) also help to stabilize the protein. Avoid repeated freezing and thawing.

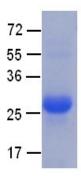

The purity of His-tagged Rab7 was determined by SDS- PAGE and Coomassie Brilliant Blue Staining.

## References

- 1. Edinger, A. L. et al., Dev. Cell 5: 571-582, 2003.
- Houlden, H. et al., Ann. Neurol. 56: 586-590, 2004.
- 3. Meggouh, F. et al., Neurology 67: 1476-1478, 2006.
- 4. Rak, A. et al., Cell 117: 749-760, 2004.
- 5. Verhoeven, K. et al., Am. J. Hum. Genet. 72: 722-727, 2003.

Tel: (610) 945-2007 Fax: (610) 945-2008

Web: www.neweastbio.com