



Product Description

— A helping hand for your research

Ras (A146P)

Catalog Number: 26478

Gene Symbol: H-Ras; K-Ras; N-Ras

Description: Anti-RAS (A146P) Mouse Monoclonal Antibody

Background: The A146P mutation of KRas results in an amino acid substitution at position 146, from an alanine to a threonine. KRAS encodes a protein that is a member of the small GTPase superfamily. Ras A146P mutation results in decreased GTPase activity and constitutive signaling. It can be found in many tumors, such as lung adenocarcinoma, mucinous adenoma, ductal carcinoma of the pancreas and colorectal carcinoma.

Immunogen: A synthetic peptide from the internal region of Ras which includes the mutation of A146P, human origin.

Tested applications: ELISA, IF, IHC

Recommended dilutions:

ELISA: 1:1000-1:5000

IF: 1:50-1:100

IHC: 1:50-1:100

Concentration: 1 mg/ml

Host: Mouse

Clonality: Monoclonal

Purity: Purified from ascites

Format: Liquid

Storage buffer:

Preservative: no

Constituents: PBS (without Mg²⁺ and Ca²⁺), pH 7.4, 150 mM NaCl, 50% glycerol


Species Reactivity: Recognizes A146P mutant, but not wild-type RAS of vertebrates.

Storage Conditions: Store at -20°C. Avoid freeze / thaw cycles.

For research use only. Not for diagnostic or therapeutic applications.

24 Whitewoods Lane, Malvern, PA 19355 Tel: (610) 945-2007 FAX: (610) 945-2008
E-mail: sale@neweastbio.com Web: www.neweastbio.com

Immunofluorescence:

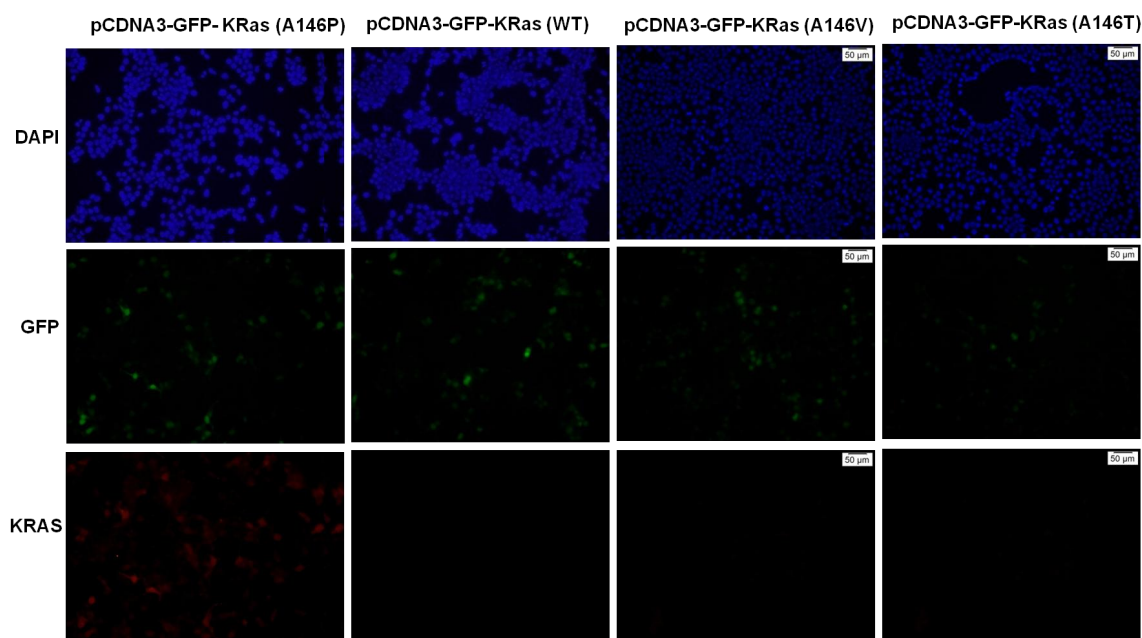

Immunofluorescence of cells expressing KRas proteins with anti-KRas (A146P) antibody.

HEK293T cells were transfected with pCDNA3-GFP-KRas (A146P) plasmid, pCDNA3-GFP-KRas (WT) plasmid, pCDNA3-GFP-KRas (A146V) plasmid or pCDNA3-GFP-KRas (A146T) plasmid, then fixed and stained with anti-KRas (A146P) monoclonal antibody (Cat. #26478).

For research use only. Not for diagnostic or therapeutic applications.

24 Whitewoods Lane, Malvern, PA 19355 Tel: (610) 945-2007 FAX: (610) 945-2008
E-mail: sale@neweastbio.com Web: www.neweastbio.com