Wuhan NewEast Biosciences Co.Ltd



Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

Arf1

Catalog Number: 26019 Gene Symbol: ARF1

Description: Anti-Arf1 Mouse Monoclonal Antibody

Background: ARF1 (ADP-ribosylation factor 1) is a member of the ARF super-family. ARF genes encode small GTPases that increase the ADP-ribosyltransferase activity of cholera toxin and are critical for vesicular trafficking as activators of phospholipase D. ARF1 protein is localized to the Golgi apparatus and has a central role in intra-Golgi transport.

Immunogen: Recombinant full length protein of

Arf1.

Tested applications: ELISA, WB, IHC, IP

Recommended dilutions:

ELISA 1:100-1:5000
WB 1:100-1:2000
IHC 1:50-1:100
Concentration: 1 mg/ml

Host: mouse

Clonality: Monoclonal

Isotype: IgG

Purity: Purified from ascites

Format: Liquid
Storage buffer:
Preservative: no


Constituents: PBS (without Mg²⁺ and Ca²⁺), pH

7.4, 150 mM NaCl, 50% glycerol

Species Reactivity: Arf1 antibody recognizes Arf1 (maybe Arf2 and Arf3) of vertebrates. **Storage Conditions:** Store at -20°C. Avoid

freeze / thaw cycles

Western blot:

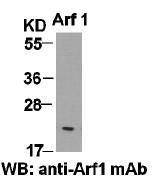

Western blot analysis of Arf1 proteins.