Wuhan NewEast Biosciences Co.Ltd



Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

Arl3

Catalog Number: 26070 Gene Symbol: ARL3

Description: Anti-Arl3 Mouse Monoclonal Antibody

Background: Arl3 (Arf-like 3) is an Arf family protein. Arl3 binds microtubules to alter specific aspects of cytokinesis. It has been proposed that retinitis pigmentosa 2 functions in concert with Arl3 to link the cell membrane and the cytoskeleton in photoreceptors as part of the cell signaling or vesicular transport machinery.

Immunogen: Recombinant Arl3 protein
Tested applications: ELISA, WB, IHC, IP

Recommended dilutions:

ELISA 1:100-1:5000 WB 1:100-1:2000 IHC 1:10-100 Concentration: 1 mg/ml

Host: Mouse

Clonality: Monoclonal

Purity: Purified from ascties

Format: Liquid Storage buffer: Preservative: no

Constituents: PBS (without Mg²⁺ and Ca²⁺), pH

7.4, 150 mM NaCl, 50% glycerol


Species Reactivity: Arl3 antibody recognizes

Arl3 from vertebrates.

Storage Conditions: Store at -20°C. Avoid

freeze / thaw cycles

Western blot:

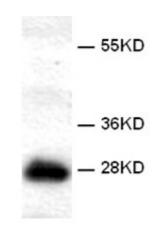

WB: Anti-Arl3 mAb

Western blot analysis of recombinant Arl3 proteins. Purified his-tagged Arl3 protein (Cat. #10152) was blotted with anti-Arl3 monoclonal antibody (Cat # 26070).