Wuhan NewEast Biosciences Co.Ltd



Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

Ga_{14}

Catalog Number: 26027 Gene Symbol: GNA14

Description: Anti-Gα₁₄ Mouse Monoclonal Antibody

Background: Heterotrimeric G proteins are essential cellular signal transducers. $G\alpha_{14}$ is a member of the Gq subfamily of G proteins. $G\alpha_{14}$ is known to stimulate phospholipase C. Other biochemical and physiological functions of

 $G\alpha_{14}$ are being explored.

Immunogen: A synthetic peptide from the

C-terminal region of human $G\alpha_{14}$.

Tested applications: ELISA, WB, IHC, IP

Recommended dilutions:

ELISA 1:100-1:5000
WB 1:100-1:2000
IHC 1:50-1:100
Concentration: 1 mg/ml

Host: mouse

Clonality: Monoclonal

Isotype: IgG

Purity: Purified from ascites

Format: Liquid Storage buffer: Preservative: no

Constituents: PBS (without Mg²⁺ and Ca²⁺), pH

7.4, 150 mM NaCl, 50% glycerol


Species Reactivity: $G\alpha_{14}$ antibody recognizes

 $G\alpha_{14}$ of vertebrates.

Storage Conditions: Store at -20°C. Avoid

freeze / thaw cycles

Western blot:

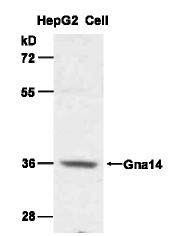

WB: anti-Gna14 mAb

 $HepG2\ whole-cell\ lysate\ was\ used\ for\ western$

blot (cat# 26027).