## Wuhan NewEast Biosciences Co.Ltd



Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

## Ran

Catalog Number: 26020 Gene Symbol: RAN

**Description:** Anti-Ran Mouse Monoclonal Antibody

**Background:** Ran (Ras-related nuclear protein) is a member of the Ras-superfamily GTPases. Ran is invovled in the transport of proteins across the nuclear envelope, as well as in microtubule organization during mitosis.

Immunogen: Recombinant full length protein of

Ran

Tested applications: ELISA, WB, IP, IHC

**Recommended dilutions:** ELISA: 1:1000-1:5000

WB: 1:100-1:2000

IP: 1 μg for 1~2 mg total cellular proteins

IHC: 1:100 dilution Concentration: 1 mg/ml

Host: mouse

Clonality: Monoclonal

Isotype: IgG

Purity: Purified from ascites

Format: Liquid Storage buffer: Preservative: no

Constituents: PBS (without Mg<sup>2+</sup> and Ca<sup>2+</sup>), pH

 $7.4,\,150$  mM NaCl, 50% glycerol


Species Reactivity: Anti-Ran antibody

recognizes Ran of vertebrates.

Storage Conditions: Store at -20°C. Avoid

freeze / thaw cycles

## Western blot:

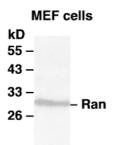

## WB: anti-Ran mAb

Western blot analysis of Ran proteins. Whole-cell lysate of MEF cells was used for western blots with the monoclonal anti-Ran antibody (Cat # 26020).