## Wuhan NewEast Biosciences Co.Ltd



Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

## BMP3

Catalog Number: 26051 Gene Symbol: BMP3

**Description:** Anti-BMP3 Mouse Monoclonal Antibody

**Background:** BMP3 belongs to the TGF-beta family. It negatively regulates bone density and antagonizes the ability of certain osteogenic BMPs to induce osteoprogenitor differentiation and ossification. BMP3 is secreted and highly expressed in fracture tissue, particularly in osteoblasts, osteoclasts and chondroblasts.

Immunogen: A synthetic peptide near the

N-terminus of BMP3, human origin.

Tested applications: ELISA, WB, IHC, IP

**Recommended dilutions:** 

ELISA: 1:1000-1:5000 WB: 1:1000-1:5000 IHC: 1:50-1:100 Concentration: 1 mg/ml

Host: Mouse

Clonality: Monoclonal

**Purity:** Purified from ascites

Format: Liquid Storage buffer: Preservative: no

Constituents: PBS (without Mg<sup>2+</sup> and Ca<sup>2+</sup>), pH

7.4, 150 mM NaCl, 50% glycerol


Species Reactivity: Anti-BMP3 antibody

recognizes BMP3 of vertebrates.

Storage Conditions: Store at -20°C. Avoid

freeze / thaw cycles.

## Western blot:

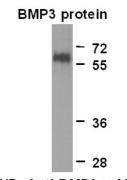

WB: Anti-BMP3 mAb

Western blot analysis of recombinant BMP3 protein. Purified His-tagged BMP3 protein (Cat. #10194) was blotted with anti-BMP3 monoclonal

antibody (Cat. #26051).