


Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

Cystatin C

Catalog Number: 26085 Gene Symbol: Cystatin C Description: Anti-Cystatin C Mouse Monoclonal Antibody

Background: Cystatin C is mainly used as a biomarker of kidney function. It is a potent inhibitor of cysteine proteases. Immunogen: Recombinant full-lemgth Cystatin C protein, human origin. Tested applications: ELISA, WB, IHC, IP **Recommended dilutions:** ELISA: 1:100-1:5000 WB: 1:100-1:1000 1:100-1:1000 IHC: Concentration: 1 mg/ml Host: Mouse **Clonality:** Monoclonal Isotype: IgG Purity: Purified from ascites Format: Liquid **Storage buffer:** Preservative: no Constituents: PBS (without Mg^{2+} and Ca^{2+}), pH 7.4, 150 mM NaCl, 50% glycerol. Species Reactivity: Cystatin C antibody recognizes Cystatin C proteins of vertebrates. Storage Conditions: Store at -20°C. Avoid freeze / thaw cycles.

Western blot:

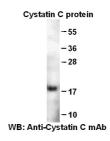

Western blot analysis of recombinant full-length Cystatin C protein. Purified His-tagged Cystatin C protein (Cat. #10157) was blotted with anti-Cystatin C monoclonal antibody (Cat. # 26085) at 1:1000 dilution.