Wuhan NewEast Biosciences Co.Ltd



Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

BRAF (V600M)

Catalog Number: 26462

Gene Symbol: BRAF, BRAF1, RAFB1

Description: Anti-BRAF (V600M) Mouse Monoclonal Antibody

Background: BRAF belongs to a family of serine-threonine protein kinases. As a part of a signaling pathway known as the RAS/MAPK pathway, it plays an important part in many cellular processes. The processes include cell proliferation, differentiation and transcriptional regulation. Mutations in the BRAF gene cause diseases. Inherited mutations in BRAF cause cardiofaciocutaneous syndrome. Acquired mutations in BRAF have been found in cancers.

Immunogen: A synthetic peptide from the internal region of BRAF which includes the mutation of V600M, human origin.

Tested applications: ELISA, IF, IHC

Recommended dilutions:

ELISA: 1:1000-1:5000 IF: 1:50-1:100 IHC: 1:50-1:100

Concentration: 1 mg/ml

Host: Mouse

Clonality: Monoclonal

Purity: Purified from ascites

Format: Liquid

Storage buffer: Preservative: no

Constituents: PBS (without Mg²⁺ and Ca²⁺), pH 7.4, 150 mM NaCl, 50% glycerol

Species Reactivity: Recognizes V600M mutant, but not wild-type BRAF of vertebrates.


Storage Conditions: Store at -20°C. Avoid freeze / thaw cycles.

Wuhan NewEast Biosciences Co.Ltd



Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

Immunofluorescence:

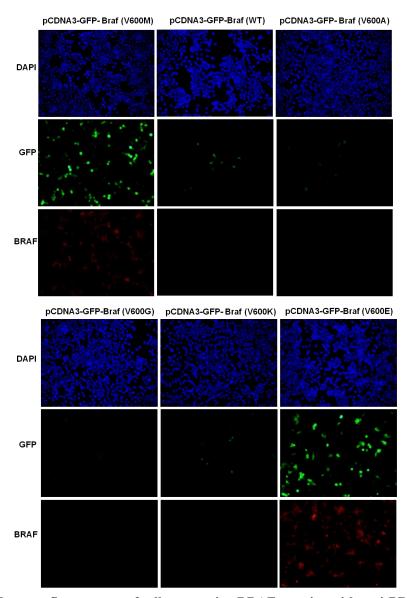

Immunofluorescence of cells expressing BRAF proteins with anti-BRAF (V600M) antibody.

HEK293T cells were transfected with pCDNA3-GFP-BRAF (V600M) plasmid, pCDNA3-GFP-BRAF (WT) plasmid,

pCDNA3-GFP-BRAF (V600A) plasmid,

pCDNA3-GFP-BRAF (V600G) plasmid,

pCDNA3-GFP-BRAF (V600K) plasmid

or pCDNA3-GFP-BRAF (V600E) plasmid

, then fixed and stained with anti-BRAF (V600M) monoclonal antibody (Cat. #26462).