## Wuhan NewEast Biosciences Co.Ltd



Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

## **TP53(R282W)**

Catalog Number: 21165 Gene Symbol: TP53

**Description:** Anti-TP53(R282W) Rabbit Polyclonal Antibody

**Background:** p53 protein acts as a tumor suppressor in many tumor types and induces growth arrest or apoptosis depending on the physiological circumstances and cell type. Defects in TP53 are found in many types of cancers.

**Immunogen:** A synthetic peptide near the C-terminus of human TP53 which includes the mutation of R282W.

Tested applications: ELISA, WB, IHC, IF

## **Recommended dilutions:**

ELISA: 1:1000-1:5000
WB: 1:100-1:1000
IHC: 1:100-1:1000
Concentration: 0.1 mg/ml

Host: Rabbit

**Clonality:** Polyclonal

Purity: Purified from serum

Format: Liquid Storage buffer: Preservative: no

Constituents: PBS (without Mg<sup>2+</sup> and Ca<sup>2+</sup>), pH

7.4, 150 mM NaCl, 50% glycerol


**Species Reactivity:** Anti-TP53(R282W) antibody recognizes TP53(R282W) of

vertebrates.

Storage Conditions: Store at -20°C. Avoid

freeze / thaw cycles

## Western blot:

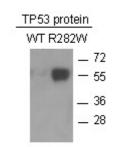

WB: Anti-TP53(R282W) pAb

Western blot analysis of recombinant TP53(R282W) and wildtype proteins. Purified His-tagged TP53(R282W) protein( Cat.#10204) and wildtype protein (Cat.#10202) were blotted with anti-TP53(R282W) rabbit antibody (Cat.#

21165).