

Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

PIK3CA (Q546R)

Catalog Number: 26214


Gene Symbol: PIK3CA, CLOVE, CWS5, MCAP, MCM, MCMTC, PI3K, p110-alpha **Description:** Anti-PIK3CA (Q546R) Mouse Monoclonal Antibody

Background: Phosphatidyl 3-kinases (PI3K) are a family of lipid kinases involved in many cellular processes, including cell growth, proliferation, differentiation, motility, and survival. It is composed of an 85 kDa regulatory subunit and a 110 kDa catalytic subunit. Mutants of PIK3CA has been found in the pathogenesis of several cancers, including colon cancer, gliomas, gastric cancer, breast cancer, endometrial cancer, and lung cancer.

Immunogen: A synthetic peptide from the internal region of PIK3CA which includes the mutation of Q546R, human origin.

Tested applications: WB Recommended dilutions: ELISA: 1:1000-1:5000 WB: 1:500-1:2000 Concentration: 1 mg/ml Host: Mouse Clonality: Monoclonal Purity: Purified from ascites Format: Liquid **Storage buffer:** Preservative: no Constituents: PBS (without Mg²⁺ and Ca²⁺), pH 7.4, 150 mM NaCl, 50% glycerol Species Reactivity: Recognizes Q546R mutant, but not wild-type PIK3CA of vertebrates. Storage Conditions: Store at -20°C. Avoid freeze / thaw cycles.

Western blot:

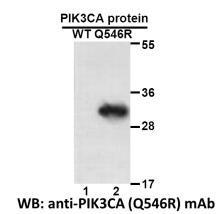

Western blot analysis of recombinant PIK3CA (Q546R) and wildtype proteins. Purified His-tagged PIK3CA (Q546R) protein (lane 2) and corresponding wildtype protein (lane 1) were blotted with anti-PIK3CA (Q546R) monoclonal antibody (Cat. #26214).