Wuhan NewEast Biosciences Co.Ltd



Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

FGFR1

Catalog Number: 26294

Gene Symbol: FGFR1, BFGFR, CD331, CEK, FGFBR, FLG **Description:** Anti-FGFR1 Mouse Monoclonal Antibody

Background: The protein FGFR1 is highly conserved between members and throughout evolution. This particular family member of FGFR1 binds both acidic and basic fibroblast growth factors and is involved in limb induction. Mutations of FGFR1 protein have been associated with Pfeiffer syndrome, Jackson-Weiss syndrome, Antley-Bixler syndrome, osteoglophonic dysplasia, squamous cell lung cancer and autosomal dominant Kallmann syndrome.

Immunogen: A synthetic peptide from the internal region of FGFR1 protein, human origin.

Tested applications: ELISA, WB, IHC

Recommended dilutions:

ELISA: 1:1000-1:2000 WB: 1:500-1:1000 IHC: 1:50-1:100 Concentration: 1 mg/ml

Host: Mouse

Clonality: Monoclonal

Purity: Purified from ascites

Format: Liquid Storage buffer: Preservative: no

Constituents: PBS (without Mg²⁺ and Ca²⁺), pH

7.4, 150 mM NaCl, 50% glycerol


Species Reactivity: Recognizes FGFR1 of

vertebrates.

Storage Conditions: Store at -20°C. Avoid

freeze / thaw cycles.

Western blot:

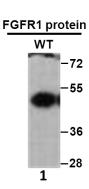

WB: Anti-FGER1 WT mAb


Western blot analysis of recombinant FGFR1 protein. Purified full length His-tagged FGFR1 (22-376aa, 45KD) protein was blotted with anti-FGFR1 monoclonal antibody (Cat. #26294).





Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

Immunofluorescence:

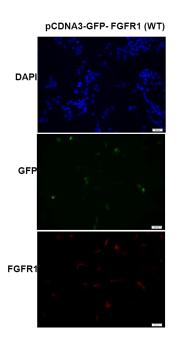

Immunofluorescence of cells expressing FGFR1 proteins with anti-FGFR1 antibody.

HEK293T cells were transfected with pCDNA3-GFP-FGFR1 (WT) (22-376) plasmid, then fixed and stained with anti-FGFR1 (WT) (22-376) monoclonal antibody (Cat. #26294).