


Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

FGFR1 (P252R)

Catalog Number: 26249 Gene Symbol: FGFR1, BFGFR, CD331, CEK, FGFBR, FLG Description: Anti-FGFR1 (P252R) Mouse Monoclonal Antibody

Background: The protein FGFR1 is highly conserved between members and throughout evolution. This particular family member of FGFR1 binds both acidic and basic fibroblast growth factors and is involved in limb induction. Mutations of FGFR1 protein have been associated with Pfeiffer syndrome, Jackson-Weiss syndrome, Antley-Bixler syndrome, osteoglophonic dysplasia, squamous cell lung cancer and autosomal dominant Kallmann syndrome. Immunogen: A synthetic peptide from the internal region of FGFR1 which includes the mutation of P252R, human origin. Tested applications: ELISA, WB, IHC **Recommended dilutions:** ELISA: 1:1000-1:2000 WB: 1:500-1:1000 IHC: 1:50-1:100 **Concentration:** 1 mg/ml Host: Mouse **Clonality:** Monoclonal **Purity:** Purified from ascites Format: Liquid Storage buffer: Preservative: no Constituents: PBS (without Mg²⁺ and Ca²⁺), pH 7.4, 150 mM NaCl, 50% glycerol Species Reactivity: Recognizes P252R mutant, but not wild-type FGFR1 of vertebrates. Storage Conditions: Store at -20°C. Avoid freeze / thaw cycles.

Western blot:

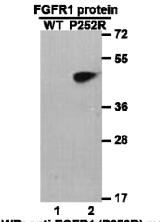

WB: anti-FGFR1 (P252R) mAb


Western blot analysis of recombinant FGFR1 (P252R) and wildtype proteins. Purified His-tagged FGFR1 (P252R) protein (22-376aa, 45KD, lane 2) and corresponding wildtype protein (lane 1) were blotted with anti-FGFR1 (P252R) monoclonal antibody (Cat. #26249).

Wuhan NewEast Biosciences Co.Ltd



Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

Immunofluorescence:

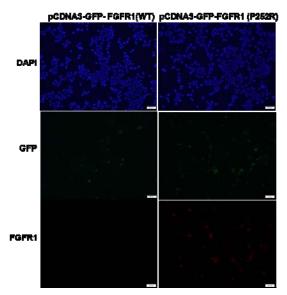

Immunofluorescence of cells expressing FGFR1 proteins with anti-FGFR1 (P252R) antibody. HEK293T cells were transfected with pCDNA3-GFP-FGFR1 (WT) plasmid (left column) or pCDNA3-GFP-FGFR1 (P252R) plasmid (right column), then fixed and stained with anti-FGFR1 (P252R) monoclonal antibody (Cat. #26249).