Wuhan NewEast Biosciences Co.Ltd



Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

RET (C634Y)

Catalog Number: 26146

Gene Symbol: RET, PTC, CDHF12, CAHF16, HSCR1, MEN2A, MEN2B

Description: Anti-RET (C634Y) Mouse Monoclonal Antibody

Background: RET protein is encoded by the RET proto- oncogene. It is a receptor tyrosine kinase. RET loss of function mutations are associated with the development of Hirschsprung's disease, while gain of function mutations are associated with development of various types of human cancer, including medullary thyroid carcinoma, multiple endocrine neoplasias type 2A and 2B, pheochromocytoma and parathyroid hyperplasia.

Immunogen: A synthetic peptide from the internal region of RET which includes the mutation of C634Y, human origin.

Tested applications: ELISA, WB, IHC

Recommended dilutions:

ELISA: 1:1000-1:2000 WB: 1:500-1:1000 IHC: 1:50-1:100

Concentration: 0.5 mg/ml

Host: Mouse

Clonality: Monoclonal **Purity:** Purified from ascites

Format: Liquid **Storage buffer:** Preservative: no

Constituents: PBS (without Mg²⁺ and Ca²⁺), pH

7.4, 150 mM NaCl, 50% glycerol


Species Reactivity: Recognizes C634Y mutant,

but not wild-type RET of vertebrates.

Storage Conditions: Store at -20°C. Avoid

freeze / thaw cycles.

Western blot:

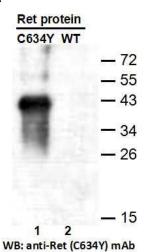

Western blot analysis of recombinant RET (C634Y) and wildtype proteins. Purified His-tagged RET (C634Y) protein (lane1) and corresponding wild-type protein (lane2) were blotted with anti-RET (C634Y) monoclonal


antibody (Cat. #26146).





Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

Immunofluorescence:

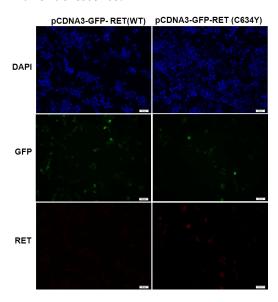

Immunofluorescence of cells expressing RET proteins with anti-RET (C634Y) antibody.

HEK293T cells were transfected with pCDNA3-GFP-RET (WT) plasmid (left column) or pCDNA3-GFP-RET (C634Y) plasmid (right column), then fixed and stained with anti-RET (C634Y) monoclonal antibody (Cat. #26146).