



Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

JAK2 (R683S)

Catalog Number: 26315

Gene Symbol: JAK2, JTK10, THCYT3

Description: Anti-JAK2 (R683S) Mouse Monoclonal Antibody

Background: JAK2 is a non-receptor tyrosine kinase. It is a number of the Janus kinase family and has been implicated in ssignaling by members of the type II cytokine receptor family, the GM-CSF receptor family, the gp130 receptor family, and the single chain receptors. JAK2 orthologs have been identified in all mammals for which complete genome genome data are available. Mutations of JAK2 have been implicated in polycythemia vera, essential thrombocythemia, and myelofibrosis as well as other myeloproliferative disorders.

Immunogen: A synthetic peptide from the internal region of JAK2 which includes the mutation of R683S, human origin.

Tested applications: ELISA, WB, IHC

Recommended dilutions:

ELISA: 1:1000-1:2000 WB: 1:500-1:1000 IHC: 1:50-1:100 Concentration: 0.5 mg/ml

Host: Mouse

Clonality: Monoclonal **Purity:** Purified from ascites

Format: Liquid **Storage buffer:** Preservative: no

Constituents: PBS (without Mg²⁺ and Ca²⁺), pH


7.4, 150 mM NaCl, 50% glycerol

Species Reactivity: Recognizes R683S mutant,

but not wild-type JAK2 of vertebrates. **Storage Conditions:** Store at -20°C. Avoid

freeze / thaw cycles.

Western blot:

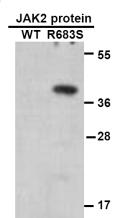

WB: anti-JAK2 (R683S) mAb


Western blot analysis of recombinant JAK2 (R683S) and wildtype proteins. Purified His-tagged JAK2 (R683S) protein (lane2) and corresponding wild-type protein (lane1) were blotted with anti-JAK2 (R683S) monoclonal antibody (Cat. #26314).

Wuhan NewEast Biosciences Co.Ltd



Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

Immunofluorescence:

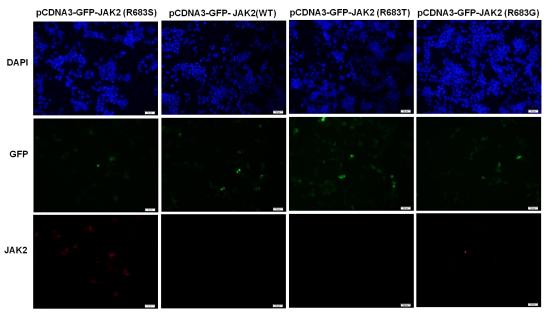

Immunofluorescence of cells expressing JAK2 proteins with anti-JAK2 (R683S) antibody. HEK293T cells were transfected with pCDNA3-GFP-JAK2 (R683S) plasmid ,pCDNA3-GFP-JAK2 (WT) plasmid, pCDNA3-GFP-JAK2 (R683T) plasmid or pCDNA3-GFP-JAK2 (R683G) plasmid, then fixed and stained with anti-JAK2 (R683S) monoclonal antibody (Cat. #26314).