#### Wuhan NewEast Biosciences Co.Ltd



Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

## KIT (V560DEL)

Catalog Number: 26289

Gene Symbol: C-Kit; CD117; KIT; PBT; SCFR

**Description:** Anti-KIT (V560DEL) Mouse Monoclonal Antibody

Background: Tyrosine-protein kinase Kit, also known as proto-oncogene c-Kit or CD117 or Mast/stem cell growth factor receptor (SCFR), is a receptor tyrosine kinase protein that in humans is encoded by the KIT gene. Activating mutations in this gene are associated with gastrointestinal stromal tumors, testicular seminoma, mast cell disease, melanoma, acute myeloid leukemia, while inactivating mutations are associated with the genetic defect piebaldism.

**Immunogen:** A synthetic peptide from the internal region of KIT which includes the deletion of V560, human origin.

Tested applications: ELISA, WB, IHC

#### **Recommended dilutions:**

ELISA: 1:1000-1:5000 WB: 1:500-1:1000 IHC: 1:50-1:100 Concentration: 1 mg/ml

**Host:** Mouse

Clonality: Monoclonal

Purity: Purified from ascites

Format: Liquid Storage buffer: Preservative: no


Constituents: PBS (without Mg<sup>2+</sup> and Ca<sup>2+</sup>), pH

7.4, 150 mM NaCl, 50% glycerol

**Species Reactivity:** Recognizes V560DEL mutant, but not wild-type KIT of vertebrates. **Storage Conditions:** Store at -20°C. Avoid

freeze / thaw cycles.

#### Western blot:

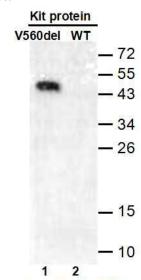

WB: anti-Kit (V560del) mAb

# Western blot analysis of recombinant KIT (V560DEL) and wild-type proteins.

Purified His-tagged KIT (V560DEL) protein (lane 1) and corresponding wild-type protein (lane 2) were blotted with anti-KIT (V560DEL)


monoclonal antibody (Cat. #26289).





Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

### Immunofluorescence:

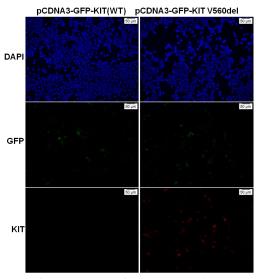

Immunofluorescence of cells expressing KIT proteins with anti-KIT (V560DEL) antibody.

HEK293T cells were transfected with pCDNA3-GFP-KIT (WT) plasmid (left column) or pCDNA3-GFP-KIT (V560DEL) plasmid (right column), then fixed and stained with anti-KIT (V560DEL) monoclonal antibody (Cat. #26289).