## Wuhan NewEast Biosciences Co.Ltd



Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

## **KIT (L576P)**

**Catalog Number: 26272** 

Gene Symbol: C-Kit; CD117; KIT; PBT; SCFR

**Description:** Anti-KIT (L576P) Mouse Monoclonal Antibody

Background: Tyrosine-protein kinase Kit, also known as proto-oncogene c-Kit or CD117 or Mast/stem cell growth factor receptor (SCFR), is a receptor tyrosine kinase protein that in humans is encoded by the KIT gene. Activating mutations in this gene are associated with gastrointestinal stromal tumors, testicular seminoma, mast cell disease, melanoma, acute myeloid leukemia, while inactivating mutations are associated with the genetic defect piebaldism.

**Immunogen:** A synthetic peptide from the internal region of KIT which includes the mutation of L576P, human origin.

Tested applications: ELISA, WB

## **Recommended dilutions:**

ELISA: 1:1000-1:5000 WB: 1:500-1:1000 Concentration: 1 mg/ml

Host: Mouse

Clonality: Monoclonal

**Purity:** Purified from ascites

Format: Liquid Storage buffer: Preservative: no

Constituents: PBS (without Mg<sup>2+</sup> and Ca<sup>2+</sup>), pH

7.4, 150 mM NaCl, 50% glycerol


Species Reactivity: Recognizes L576P mutant,

but not wild-type KIT of vertebrates.

Storage Conditions: Store at -20°C. Avoid

freeze / thaw cycles.

## Western blot:

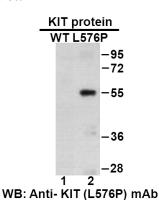

Western blot analysis of recombinant KIT

**(L576P) and wildtype proteins.** Purified His-tagged KIT (L576P) protein (lane 2) and corresponding wildtype protein (lane 1) were blotted with anti-KIT (L576P) monoclonal

antibody (Cat. #26272).