



Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

BRCA2 (2793-3030)

Catalog Number: 26364-4

Gene Symbol: BRCA2, BRCC2, FAD1,FCD, BROVCA2, PNCA2 **Description:** Anti-BRCA2 (2793-3030) Mouse Monoclonal Antibody

Background: BRCA2 is a human tumor suppressor gene. Its protein, also called by the synonym breast cancer type 2 susceptibility protein, is responsible for repairing DNA. Mutations of BRCA2 can cause an increased risk of breast cancer. In addition to breast cancer in men and women, mutations in BRCA2 also lead to an increased risk of ovarian, Fallopian tube, prostate, and pancreatic cancers, as well as malignant melanoma.

Immunogen: Purified BRCA2 (2793-3030)

protein, human origin.

Tested applications: ELISA, WB, IHC

ELISA: 1:1000-1:2000 WB: 1:500-1:1000 IHC: 1:50-1:100 Concentration: 1 mg/ml

Recommended dilutions:

Host: Mouse

Clonality: Monoclonal **Purity:** Purified from ascites

Format: Liquid
Storage buffer:
Preservative: no


Constituents: PBS (without Mg²⁺ and Ca²⁺), pH

7.4, 150 mM NaCl, 50% glycerol **Species Reactivity:** Recognizes BRCA2(2793-3030) of vertebrates.

Storage Conditions: Store at -20°C. Avoid

freeze / thaw cycles.

Western blot:

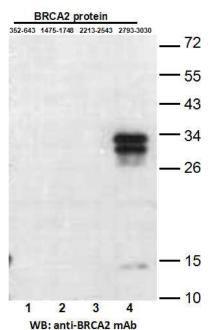


Western blot analysis of recombinant BRCA2 proteins. Purified His-tagged BRCA2 (352-643) protein (lane1), His-tagged BRCA2 (1475-1748) protein(lane2), His-tagged BRCA2 (2213-2543) protein(lane3) and His-tagged BRCA2 (2793-3030) protein (lane4) were blotted with anti-BRCA2 (2793-3030) monoclonal antibody (Cat. #26364-4).

Wuhan NewEast Biosciences Co.Ltd



Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

Immunofluorescence:

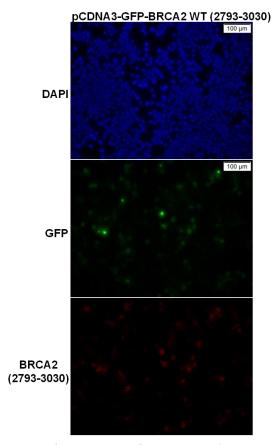

Immunofluorescence of cells expressing BRCA2 (2793-3030) protein with anti-BRCA2 (2793-3030) antibody. HEK293T cells was transfected with pCDNA3-GFP-BRCA2 (2793-3030) plasmid, then fixed and stained with anti-BRCA2 (2793-3030) monoclonal antibody (Cat. #26364-4).