



Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

ERBB2 (D769H)

Catalog Number: 26381

Gene Symbol: ERBB2/HER2/Neu/c-Neu

Description: Anti-ERBB2 (D769H) Mouse Monoclonal Antibody

Background: The HER-2/neu oncogene, a member of the epidermal growth factor receptor or erb-B gene-like family, encodes a transmembrane tyrosine kinase receptor that mediates extracellular signals activated by epidermal growth factors. Her2 abnormal has been strongly associated with many malignant tumors, especially with breast cancers. The expression level of Her2 is an important critiria in clinic evaluating of the progression of breast cancer.

Immunogen: A synthetic peptide from the internal region of ERBB2 which includes the mutation of D769H, human origin.

Tested applications: ELISA, WB, IHC

Recommended dilutions:

ELISA: 1:1000-1:2000 WB: 1:500-1:1000 IF: 1:25-1:50

Concentration: 1 mg/ml

Host: Mouse

Clonality: Monoclonal

Purity: Purified from ascites

Format: Liquid Storage buffer: Preservative: no

Constituents: PBS (without Mg²⁺ and Ca²⁺), pH


7.4, 150 mM NaCl, 50% glycerol

Species Reactivity: Recognizes D769H mutant,

but not wild-type ERBB2 of vertebrates. **Storage Conditions:** Store at -20°C. Avoid

freeze / thaw cycles.

Western blot:

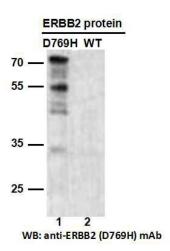

Western blot analysis of recombinant ERBB2 (D769H) and wildtype proteins. Purified His-tagged ERBB2 (D769H) protein (lane 1) and corresponding wildtype protein (lane 2) were blotted with anti-ERBB2 (D769H) monoclonal


antibody (Cat. #26381).





Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

Immunofluorescence:

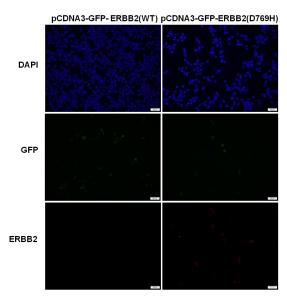

Immunofluorescence of cells expressing ERBB2 proteins with anti-ERBB2 (D769H) antibody. HEK293T cells were transfected with pCDNA3-GFP-ERBB2 (WT) plasmid (left column) or pCDNA3-GFP-ERBB2 (D769H) plasmid (right column), then fixed and stained with anti-ERBB2 (D769H) monoclonal antibody (Cat. #26381).