


Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

CTNNB1 (T41A)

Catalog Number: 26305 Gene Symbol: Beta-catenin, CTNNB, Catenin beta-1 Description: Anti-CTNNB1 (T41A) Mouse Monoclonal Antibody

Background: CTNNB1 protein is a dual function protein. It is a subunit of a complex of proteins that from adherens junctions, which are important for the establishment and maintenance of epithelial cell layers by regulating cell growth and adhesion between adjacent cells. CTNNB1 protein also pulls double duty as an intracellular signal transducer in the Wnt signaling pathway. Mutations of CTNNB1 have been implicated in the pathogenesis of several cancers. Immunogen: A synthetic peptide from the internal region of CTNNB1 which includes the mutation of T41A, human origin. Tested applications: ELISA, WB, IHC **Recommended dilutions:** ELISA: 1:1000-1:5000 WB∙ 1:500-1:2000 IHC: 1:50-1:100 **Concentration:** 1 mg/ml Host: Mouse **Clonality:** Monoclonal Purity: Purified from ascites Format: Liquid **Storage buffer:** Preservative: no Constituents: PBS (without Mg²⁺ and Ca²⁺), pH 7.4, 150 mM NaCl, 50% glycerol Species Reactivity: Recognizes T41A mutant, but not wild-type CTNNB1 of vertebrates. Storage Conditions: Store at -20°C. Avoid freeze / thaw cycles.

Western blot:

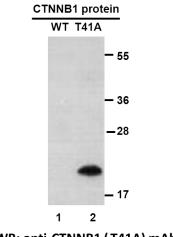

WB: anti-CTNNB1 (T41A) mAb


Western blot analysis of recombinant CTNNB1 (T41A) and wildtype proteins. Purified His-tagged CTNNB1 (T41A) protein (lane2) and corresponding wild-type protein (lane1) were blotted with anti-CTNNB1 (T41A) monoclonal antibody (Cat. #26305).

Wuhan NewEast Biosciences Co.Ltd



Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

Immunofluorescence:

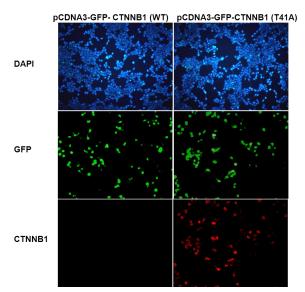

Immunofluorescence of cells expressing CTNNB1 proteins with anti-CTNNB1 (T41A) antibody. HEK293T cells were transfected with pCDNA3-GFP-CTNNB1 (WT) plasmid (left column) or pCDNA3-GFP-CTNNB1 (T41A) plasmid (right column), then fixed and stained with anti-CTNNB1 (T41A) monoclonal antibody (Cat. #26305).