Wuhan NewEast Biosciences Co.Ltd



Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

Rab3A

Catalog Number: 21041 Gene Symbol: RAB3A

Description: Anti-Rab3A Rabbit Polyclonal Antibody

Background: Rab3A is a member of the Rab-family GTPases and is specifically targeted

to synaptic vesicles.

Immunogen: A synthetic peptide near the

C-terminus of human Rab3A.

Tested applications: ELISA, WB, IHC, IP

Recommended dilutions:

ELISA 1:100-1:5000
WB 1:100-1:2000
IHC 1:50-100
Concentration: 0.2 mg/ml

Host: Rabbit

Clonality: Polyclonal

Purity: Purified from serum

Format: Liquid **Storage buffer:** Preservative: no

Constituents: PBS (without Mg²⁺ and Ca²⁺), pH

7.4, 150 mM NaCl, 50% glycerol

Species Reactivity: Rab3A antibody recognizes


Rab3A of vertebrates.

Storage Conditions: Store at -20°C. Avoid

freeze / thaw cycles

Western blot:

GST-Rab3A

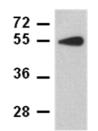

WB: anti-Rab3A pAb

Western blot analysis of anti-Rab3A antibody. Purified GST-tagged Rab3A protein was used for Western blot with anti-Rab3A polyclonal

antibody (Cat. # 21041).