



Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

RheB

Catalog Number: 21098 Gene Symbol: RHEB

Description: Anti-RheB Rabbit Polyclonal Antibody

Background: Small GTPases are a super-family of cellular signaling regulators. RheB is a member of the Ras-superfamily GTPases. RheB has been shown to interact with C-Raf, Mammalian target of rapamycin (mTOR), TSC2, Ataxia telangiectasia mutated (ATM), KIAA1303 and Ataxia telangiectasia and Rad3

related (ATR).

Immunogen: A synthetic peptide from the

C-terminal region of human RheB **Tested applications:** ELISA, WB, IHC

Recommended dilutions:

ELISA: 1:1000-1:5000 WB: 1:100-1:1000 IHC: 1:100-1:1000 Concentration: 0.2 mg/ml

Host: Rabbit

Clonality: Polyclonal

Purity: Purified from antisera

Format: Liquid Storage buffer: Preservative: no

Constituents: PBS (without Mg²⁺ and Ca²⁺), pH

7.4, 150 mM NaCl, 50% glycerol


Species Reactivity: Anti-RheB antibody

recognizes RheB of vertebrates.

Storage Conditions: Store at -20°C. Avoid

freeze / thaw cycles

Western blot:

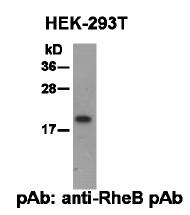

Western blot analysis of RheB proteins. Whole-cell lysate of HEK-293T cells was used for western blot with the polyclonal anti-RheB antibody (Cat # 21098).