



Email:neweastbio@126.com Web:www.neweastbio.com.cn Tel:027-87561033

Ran

Catalog Number: 21097 Gene Symbol: RAN

Description: Anti-Ran Rabbit Polyclonal Antibody

Background: Ran (Ras-related nuclear protein) is a member of the Ras-superfamily GTPases. Ran is invovled in the transport of proteins across the nuclear envelope, as well as in microtubule organization during mitosis.

Immunogen: A synthetic peptide from the

C-terminal region of human Ran

Tested applications: ELISA, WB, IHC

 Recommended dilutions:

 ELISA:
 1:1000-1:5000

 WB:
 1:100-1:2000

 IHC:
 1:100-1:1000

 Concentration:
 0.2 mg/ml

Host: Rabbit

Clonality: Polyclonal

Purity: Purified from antisera

Format: Liquid **Storage buffer:** Preservative: no

Constituents: PBS (without Mg²⁺ and Ca²⁺), pH


7.4, 150 mM NaCl, 50% glycerol **Species Reactivity:** Anti-Ran antibody

recognizes Ran of vertebrates.

Storage Conditions: Store at -20°C. Avoid

freeze / thaw cycles

Western blot:

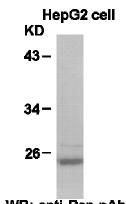

WB: anti-Ran pAb

Western blot analysis of Ran proteins. Whole-cell lysate of HepG2 cells was used for western blots with the polyclonal anti-Ran antibody (Cat # 21097).